


Up to 20% of patients will require re-intervention following iliofemoral venous stenting.¹

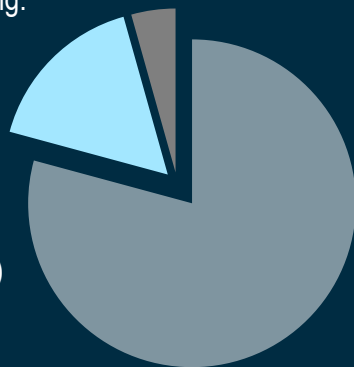

In-stent restenosis (ISR) accounts for the majority (≤ 83%) of re-interventions.¹

Crossing in-stent occlusions is not always possible with conventional techniques.²⁻⁴

POWER up your workflow.

In cases where standard techniques had FAILED

Operators **SUCCESSFULLY RECANALIZED 67- 89%** of chronically occluded venous stents with the **PowerWire®** family of RF Guidewires.²⁻⁴

1. Saleem, T. and R. Seshadri. (2022). An overview of in-stent restenosis in iliofemoral venous stents. *JVS-VL*, 10(2):P492-503.
2. Neidert, N., and H. Bjarnson. (2019). Abstract: Effectiveness of the PowerWire Radiofrequency Guidewire in Recanalizing Chronically Occluded Iliac Venous Stents. *JVS-VL*; 7(2):299-300.
3. Majdalany, B.S., et al. (2018). Radiofrequency Wire Recanalization of Chronically Occluded Venous Stents: A Retrospective, Single-Center Experience in 15 Patients. *Cardiovasc Intervent Radiol*; 42(1):130-136.
4. Shapiro, J., et al. (2022). Novel therapy for recanalization of chronic iliofemoral venous occlusion using radiofrequency. *J Vasc Curg Venous Lymphatic Disord*; 10(6):1288-1293.

Product Specifications

Cather Compatibility	4F (minimum)
Maximum OD	0.035"
Length	250 cm
Tip configurations	Straight 20°, 30°, 40° Angled

Ordering Information

Product Code	Description	Tip Shape & Strength
PSK35-250-10-6S	PowerWire® Pro RF Guidewire Kit 75 Straight	Straight, 75g
PSK35-250-12-6S	PowerWire® Pro RF Guidewire Kit 110 Straight	Straight, 110g
PSK35-250-12-6A-20	PowerWire® Pro RF Guidewire Kit Angled 20	Angled, 110g
PSK35-250-12-6A-30	PowerWire® Pro RF Guidewire Kit Angled 30	Angled, 110g
PSK35-250-12-6A-40	PowerWire® Pro RF Guidewire Kit Angled 40	Angled, 110g



Baylis
MED TECH

Baylis Medical Technologies Inc.
2645 Matheson Blvd East
Mississauga, ON Canada L4W 5S4

Tel: 1 (888)-505-4885
www.baylismedtech.com


*The PowerWire® Pro RF Guidewire is cleared by FDA to create a channel in totally occluded peripheral vessels 3 mm or greater, including vessels with stents.

PRM-00905 EN J-1 V-1 © Copyright Baylis Medical Technologies Inc., 2023. The Baylis Medical Technologies logo, Advancing Technology. Transforming Care., and PowerWire are trademarks or registered trademarks of Baylis Medical Technologies Inc. in the United States and/or other countries. Baylis Medical Technologies Inc. reserves the right to change specifications or to incorporate design changes without notice and without incurring any obligation relating to equipment previously manufactured or delivered. Patents Pending and/or issued. CAUTION: Federal Law (USA) restricts the sale of these devices to or by the order of a physician. Before use, consult product labels and inserts for any indications, contraindications, hazards, warnings, cautions and instructions for use.

Products shown may not be approved or available for sale in all jurisdictions.


Baylis Advancing Technology.
MED TECH Transforming Care.™

Challenging Occlusion?
Power through.

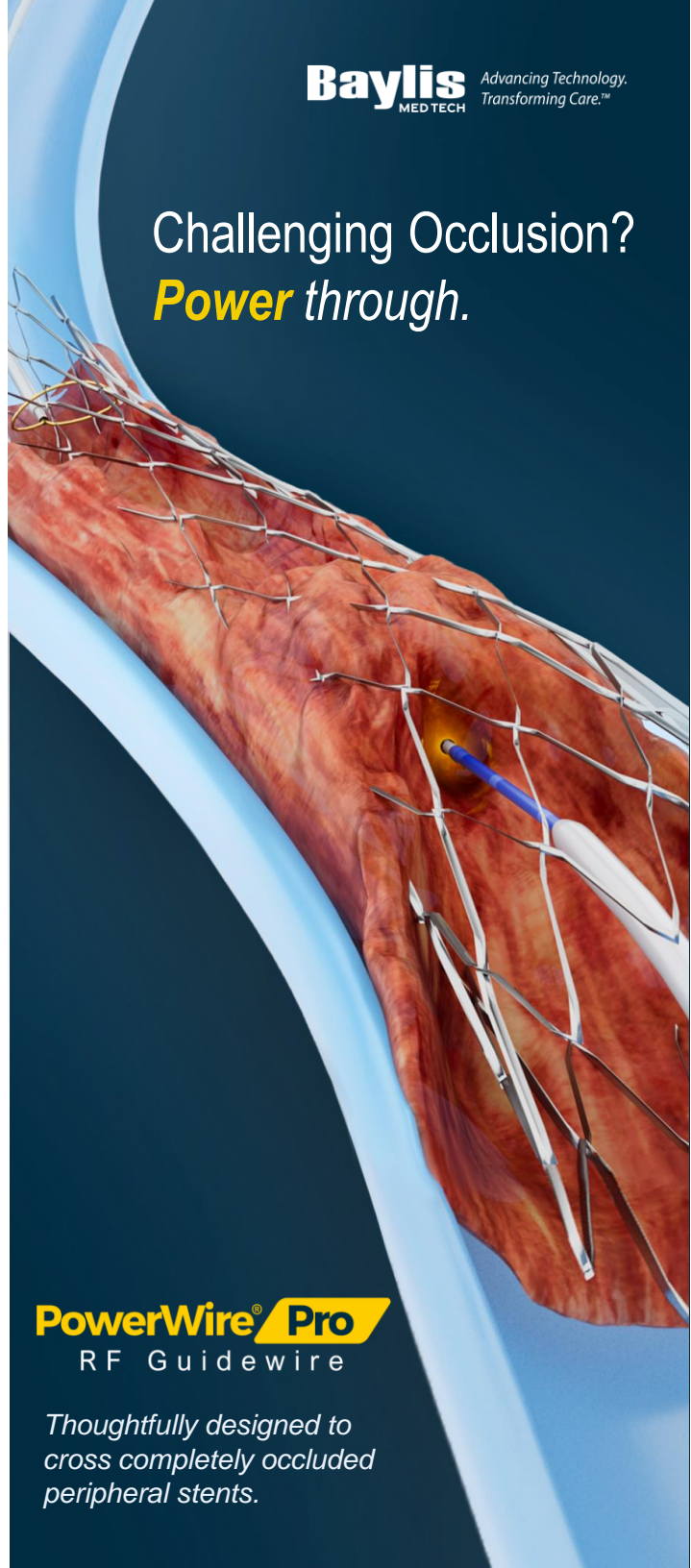

PowerWire® Pro
RF Guidewire

Thoughtfully designed to cross completely occluded peripheral stents.

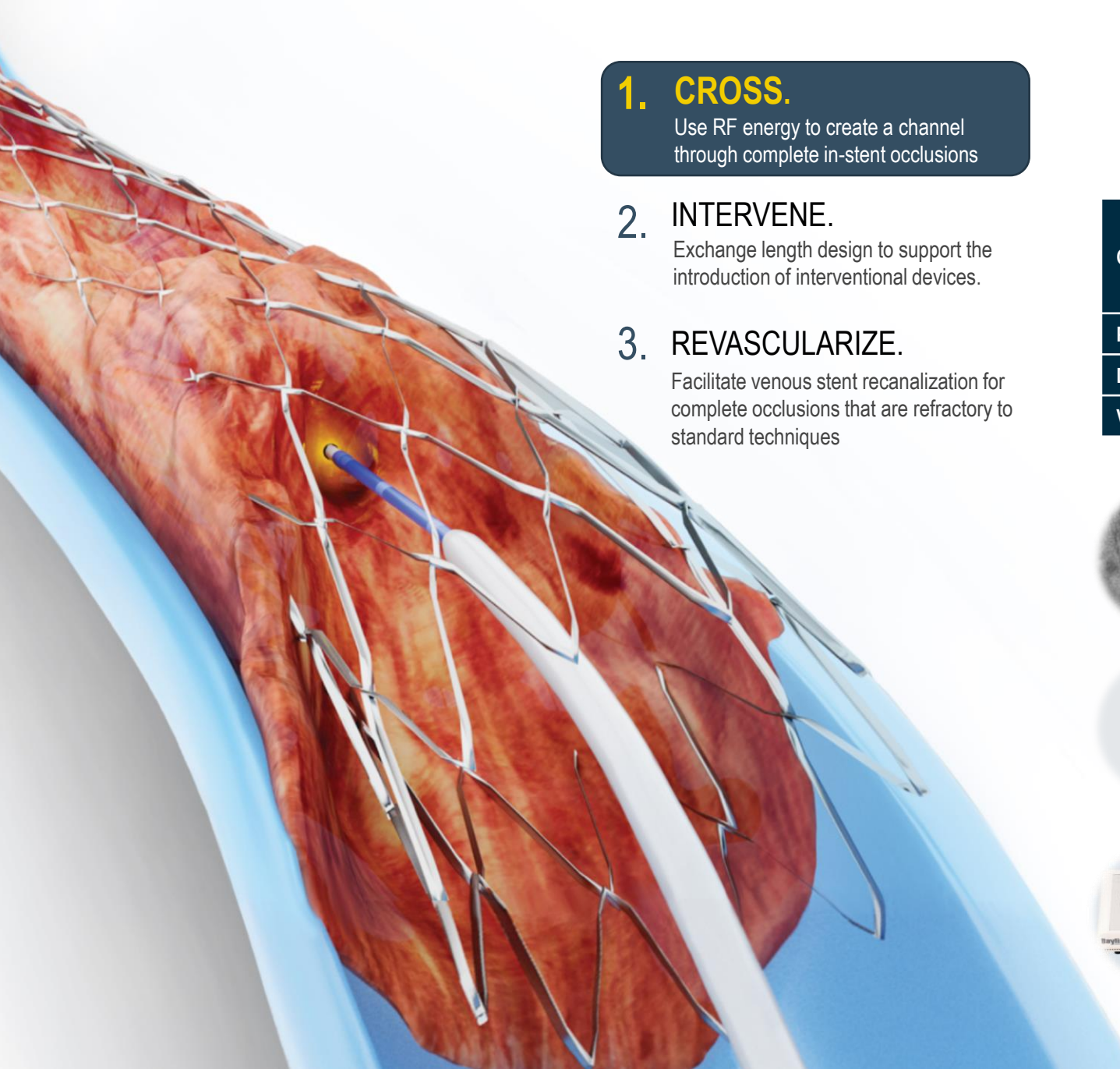

1. CROSS.

Use RF energy to create a channel through complete in-stent occlusions

2. INTERVENE.

Exchange length design to support the introduction of interventional devices.

3. REVASCULARIZE.


Facilitate venous stent recanalization for complete occlusions that are refractory to standard techniques

Cross in-stent occlusions with **RF PUNCTURE** technology

(✓) RF Puncture

(X) RF Ablation

Objective	Create a small opening with minimal damage to surrounding tissue	Create a lesion with thermal destruction of surround tissue
Power	Low (5-25W)	High (30-50W)
Duration	Short (0.3-3 s)	Long (60-90 s)
Voltage	High (270-400V)	Low (35-50V)

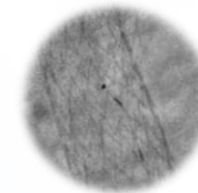

VISUALIZE the radiopaque tip within the stent under fluoroscopy



Selectively **APPLY RF ENERGY** to cross segments of the occlusion that cannot be traversed mechanically



Contact with metal terminates RF energy, **REDUCING THE POTENTIAL FOR VESSEL EXTRAVASATION** when crossing in-stent occlusions

CREATE A PATH FORWARD

Enable intervention with reliable crossing.

PowerWire[®] Pro
RF Guidewire